




**www.medsourcerental.com**

***MEDSource Inc.***

## **Wall Mounted Light**

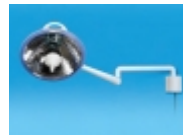

The wall mounted models are available for the MH Surgery Light, Video Monitor Arm, or Video Camera, allowing for maximum flexibility, where wall mounting is required.

The next level in the evolution of surgery lighting has been achieved with the new MH (High Intensity Discharge) Series. This all new design heralds a notable increase in efficiency over traditional halogen light sources, producing over 100,000 lux output with only a 39 watt H.I.D. bulb, versus a 150 watt halogen bulb. Heat has been virtually eliminated, while color temperature has been significantly increased to over 4,300° Kelvin. Additionally, the MH light boasts a rated bulb life of over 10,000 hours and includes complete mechanical dimming, the most accurate method.

- Most efficient surgery lighting source on the market
- Extremely cool output - lowest possible heat
- Superb color temperature
- Extremely long bulb life
- Wide range mechanical dimming

Additional Options:

Monitor Arm

- Sterile handle for adjustments within the sterile field
- Provides options for High Definition, fiber optics cabling, video source switching, and compatibility with PACS
- Future mountings of dual monitors or monitor/camera will be available
- Can be ordered with or without monitor

- Ceiling mount version as solo arm or shared with lights
- Wall mounted model also available

#### P1 Video Camera

- Provides "S-Video" imaging
- 27x optical zoom
- Auto and manual focus
- Sterilizable handle
- Mounted on its own independent arm assembly
- Can be ordered in combination with one or two lights on a Duo or Trio ceiling mounting

#### [Vendor Information](#)